




WOODPEX X

Precision Beyond Measure



APEX LOCATOR WITH PULP TESTER

New Circuit Design Minimized Interference

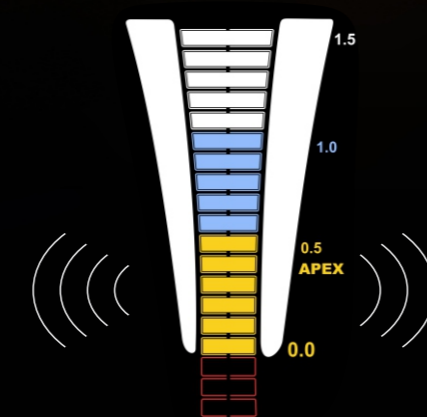

Sound Recognition Capable Aid in Operation

Two sound modes
Adjust volume and mode as you wish



Working Time: 7h Charging Time: 3h Battery: 3.7V/2000mAh Display: 3.8-inch LCD screen

Large Color Screen Clearer and Visual Display

3.8-inch LCD screen

Magnification to display the apical third

Clinical recommendation:

Measure to the 0.0 position and manually
subtract 0.5~1mm for the final working length

Four Working Modes


6.2s to reach the maximum pulse
intensity at the quickest

Mode 1: 42.5s Mode 2: 32s
Mode 3: 24s Mode 4: 6.2s

Support Apical Stop Reset

The apical stop can be freely set according
to clinical needs.

Reset the display point that is 0.5mm from the apex.
Provide warning prompts based on the
preset apical stop.

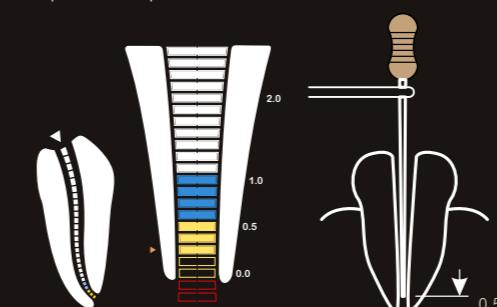

Added Sleep Mode Automatic Wake-up

Three sleep durations available

The device will be woken up upon the
measurement in sleep mode

10s 20s 30s

Comprehensive Configuration

Not autoclavable

Measuring wire x 1 (Soft texture)


Autoclavable

File clip x 4 (Coated with biological silica gel)

Lip hook x 5 (Made of medical-grade stainless steel)

Apex locator probe x 2

Pulp tester probe x 2

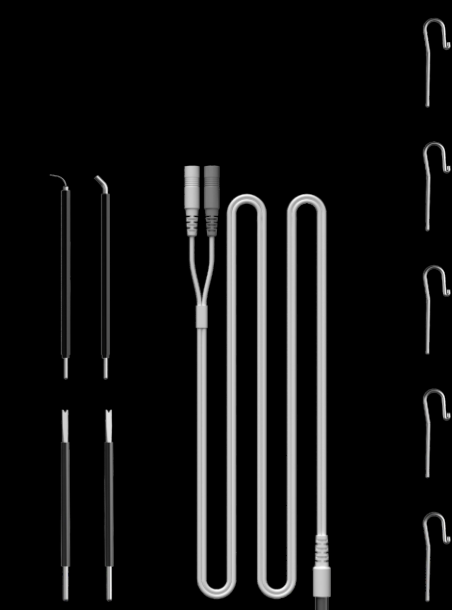

Guilin Woodpecker Medical Instrument Co., Ltd.
Information Industrial Park, Guilin National High-Tech
Zone, Guilin, Guangxi, 541004 P.R.China

Tel: +86-773-5873196
Fax: +86-773-5822450
E-mail: woodpecker4@glwoodpecker.com
Website: www.glwoodpecker.com

ECREP MedNet GmbH
Borkstrasse 10 · 48163 Muenster · Germany

DISTRIBUTED BY:

